

MKM Transparent flask

Developed by Dipl.zt. Martin Mikuš



Transparent MKM flask

Enhancements:

Uniform polymerization

Transparent for gingival toning

PMMA does not adhere to the flask


Sized for smaller polymerization units

Saves duplicating material


Also for making stone and investment duplicate models



Comparison of castings

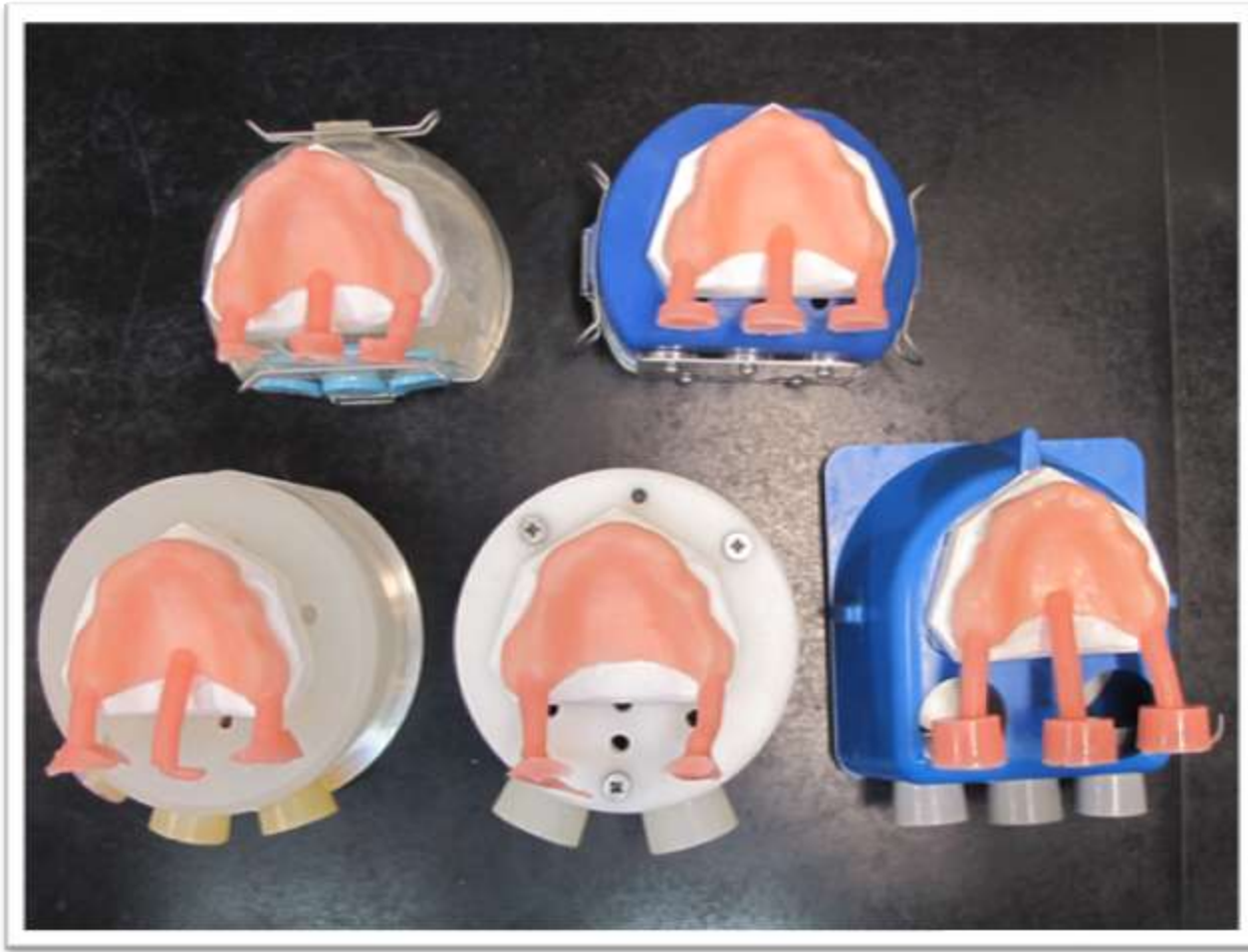

MKM Transparent

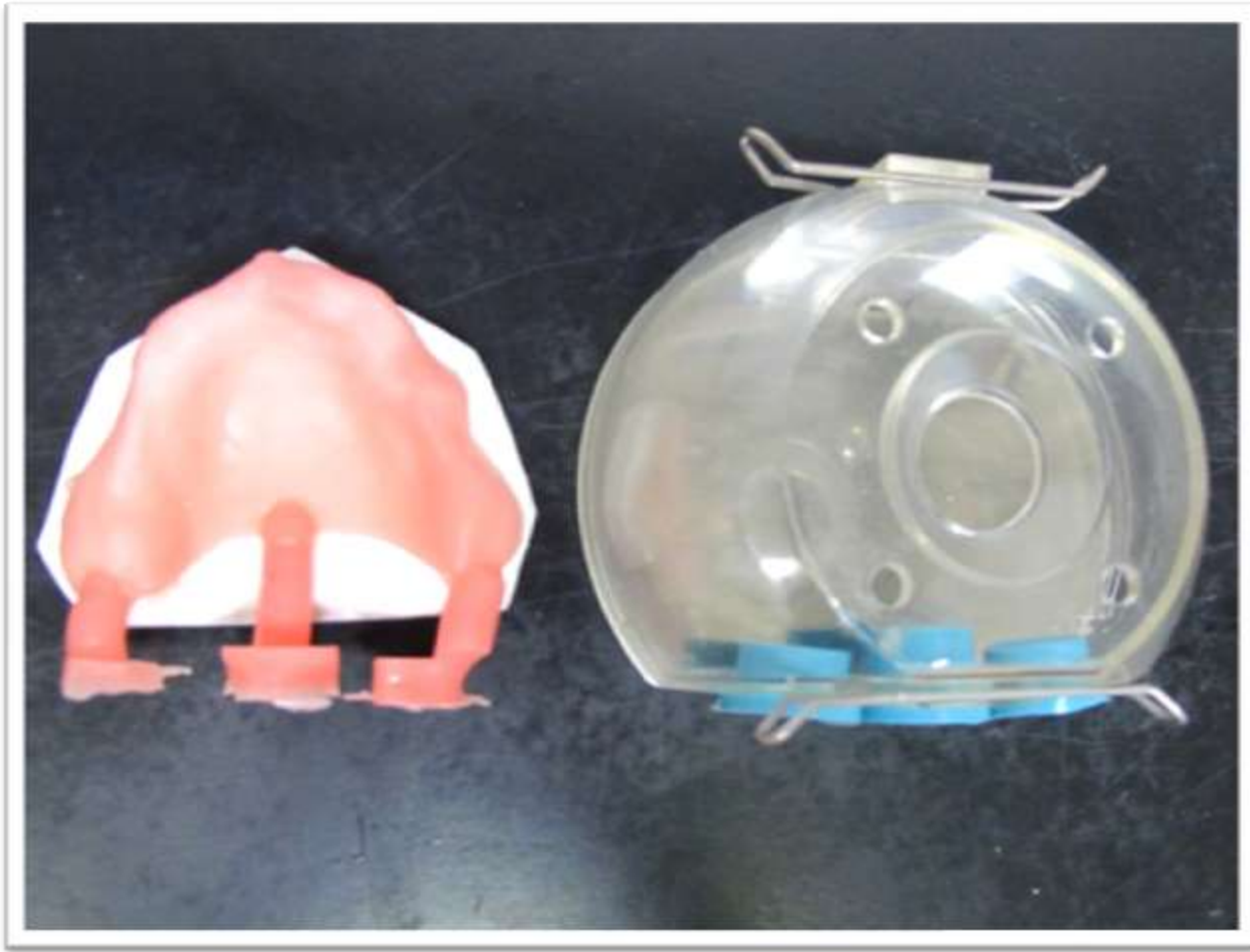


MKM Transparent




MKM Transparent




Type 1

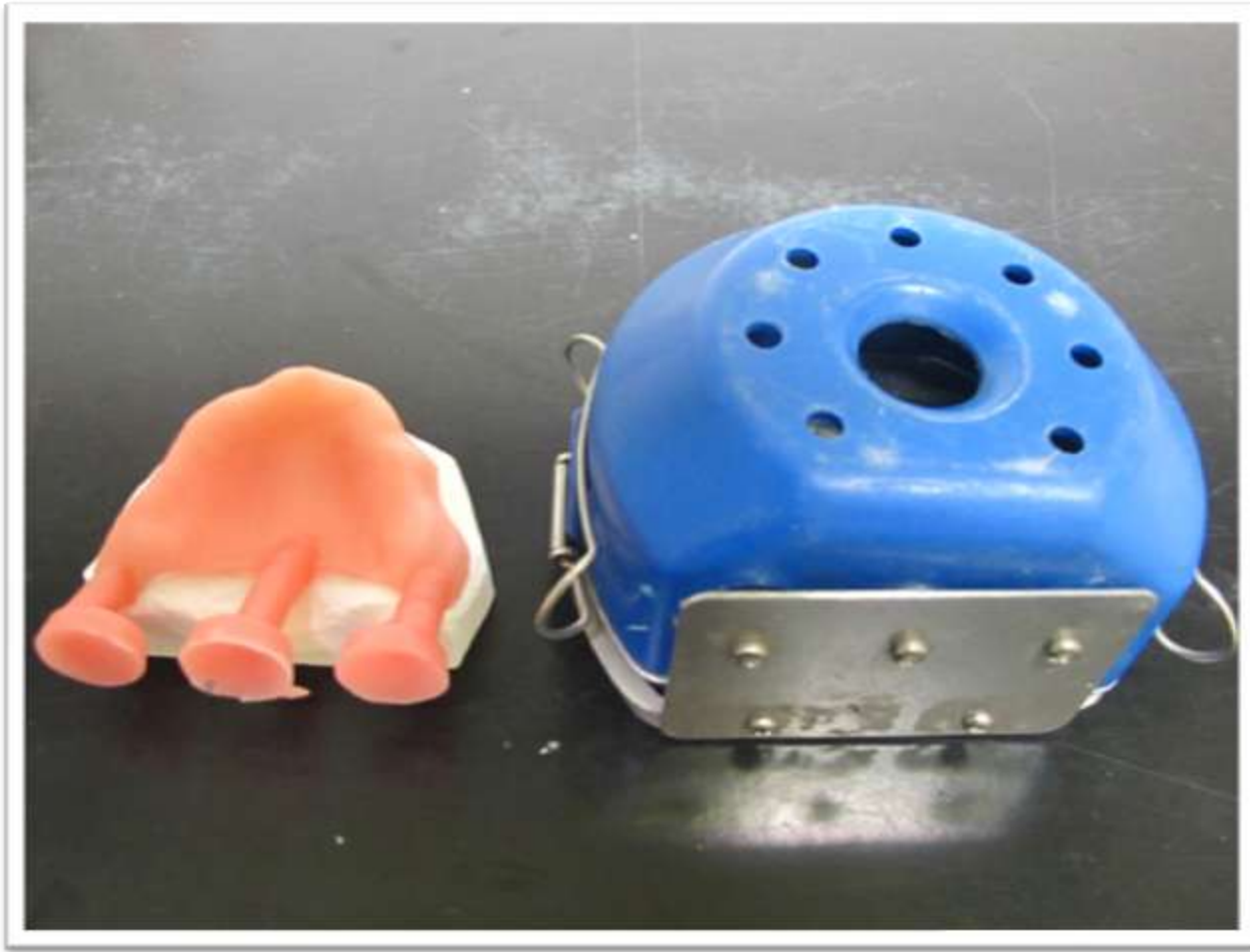


Type 1

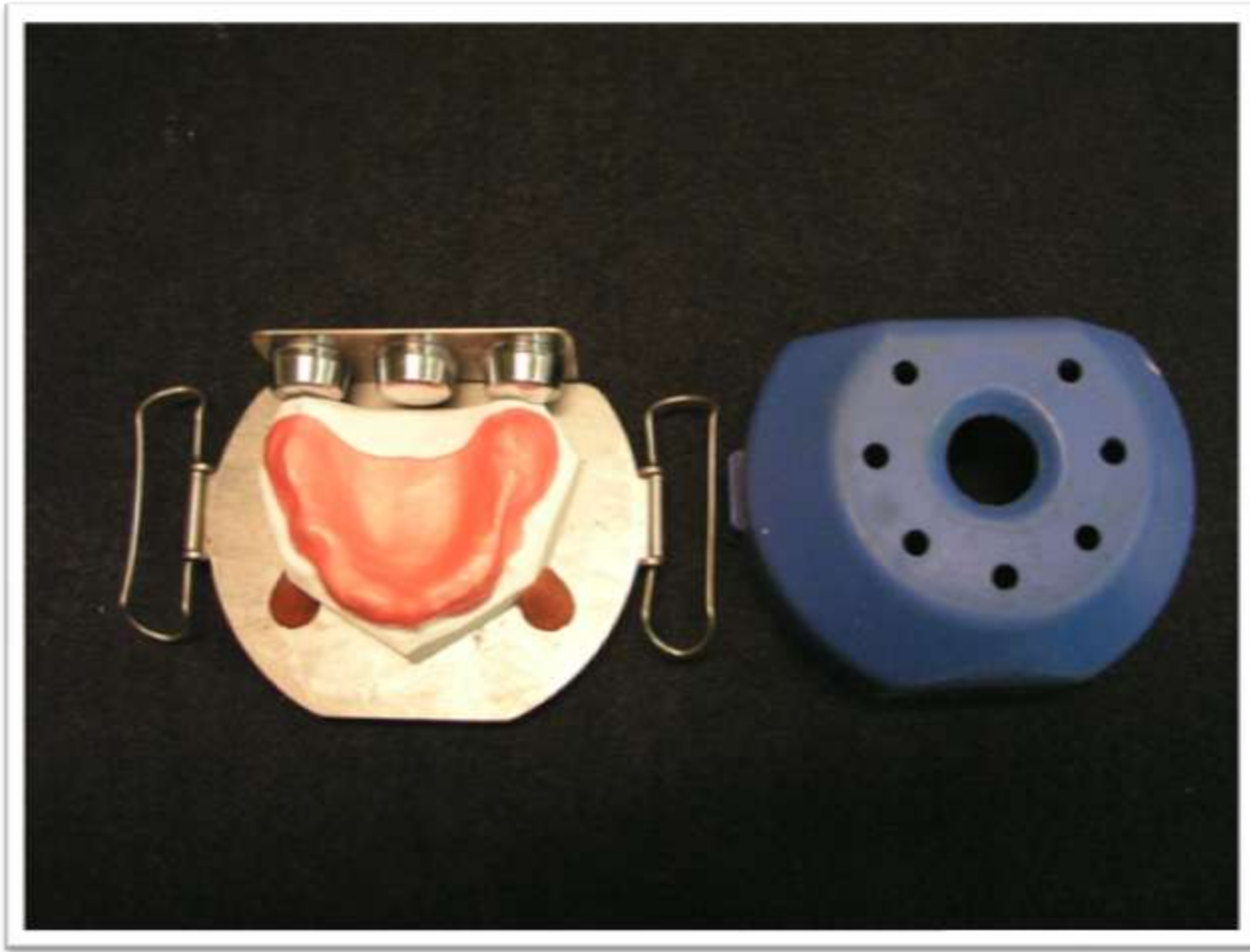


Type 1

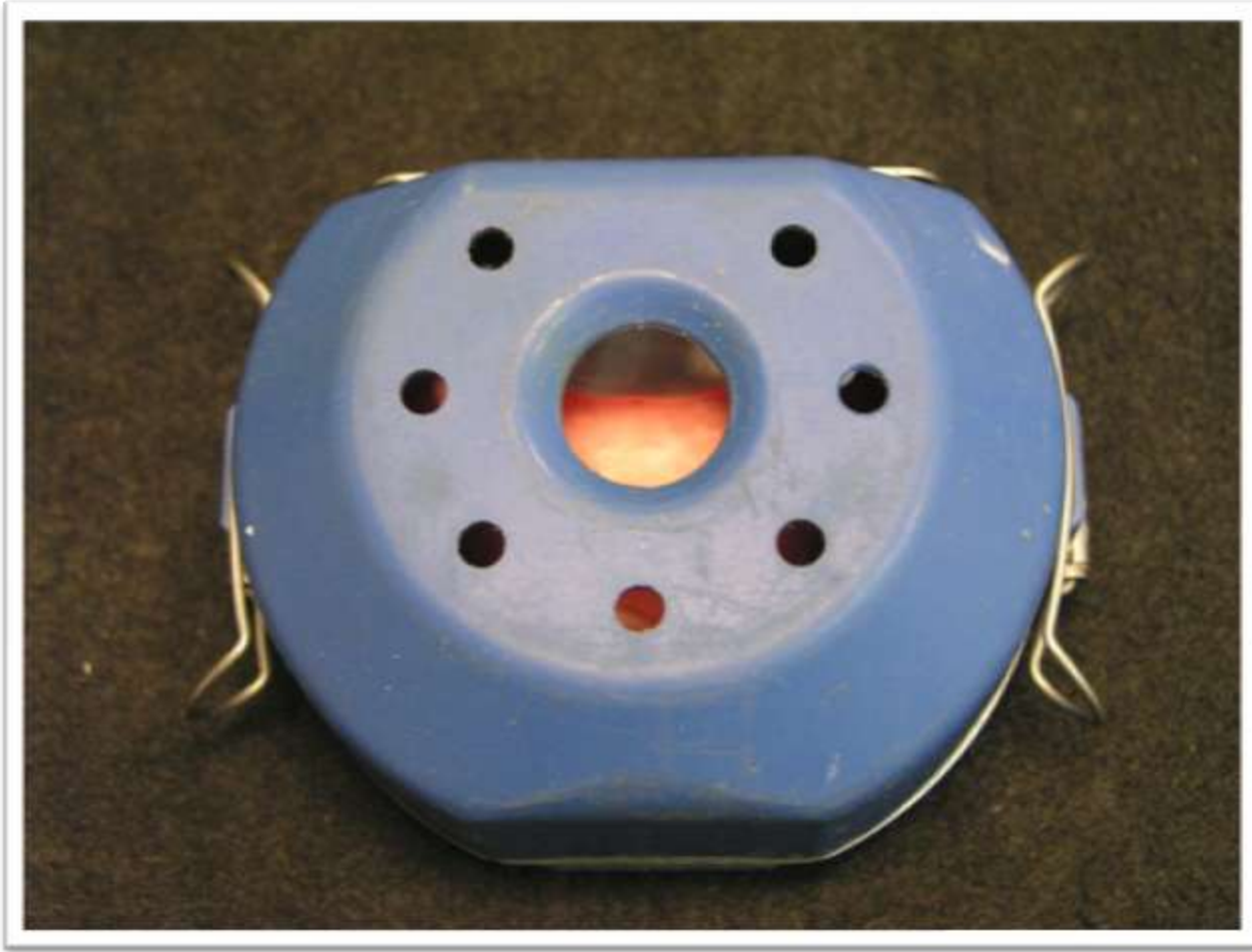

Type 2

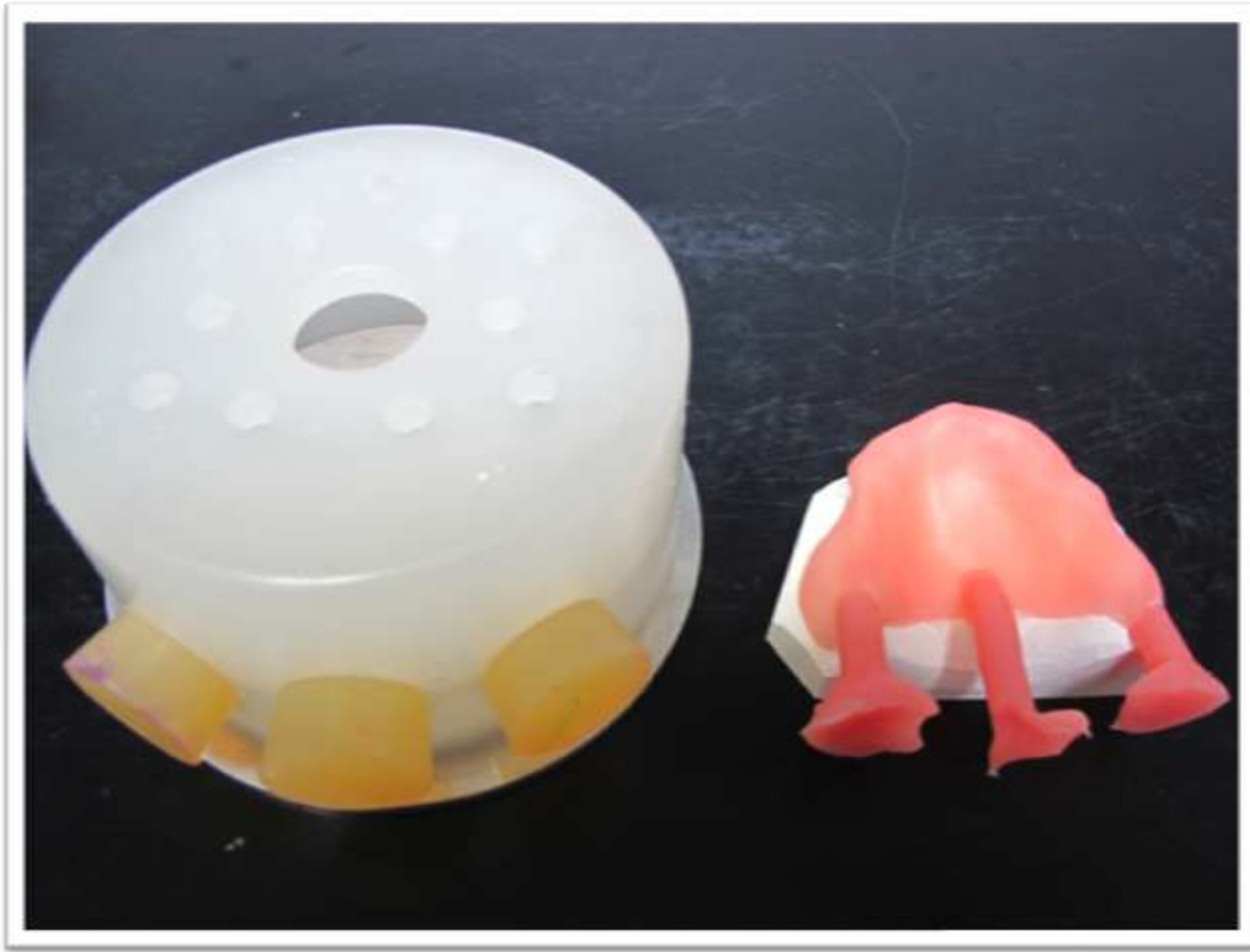


Type 2

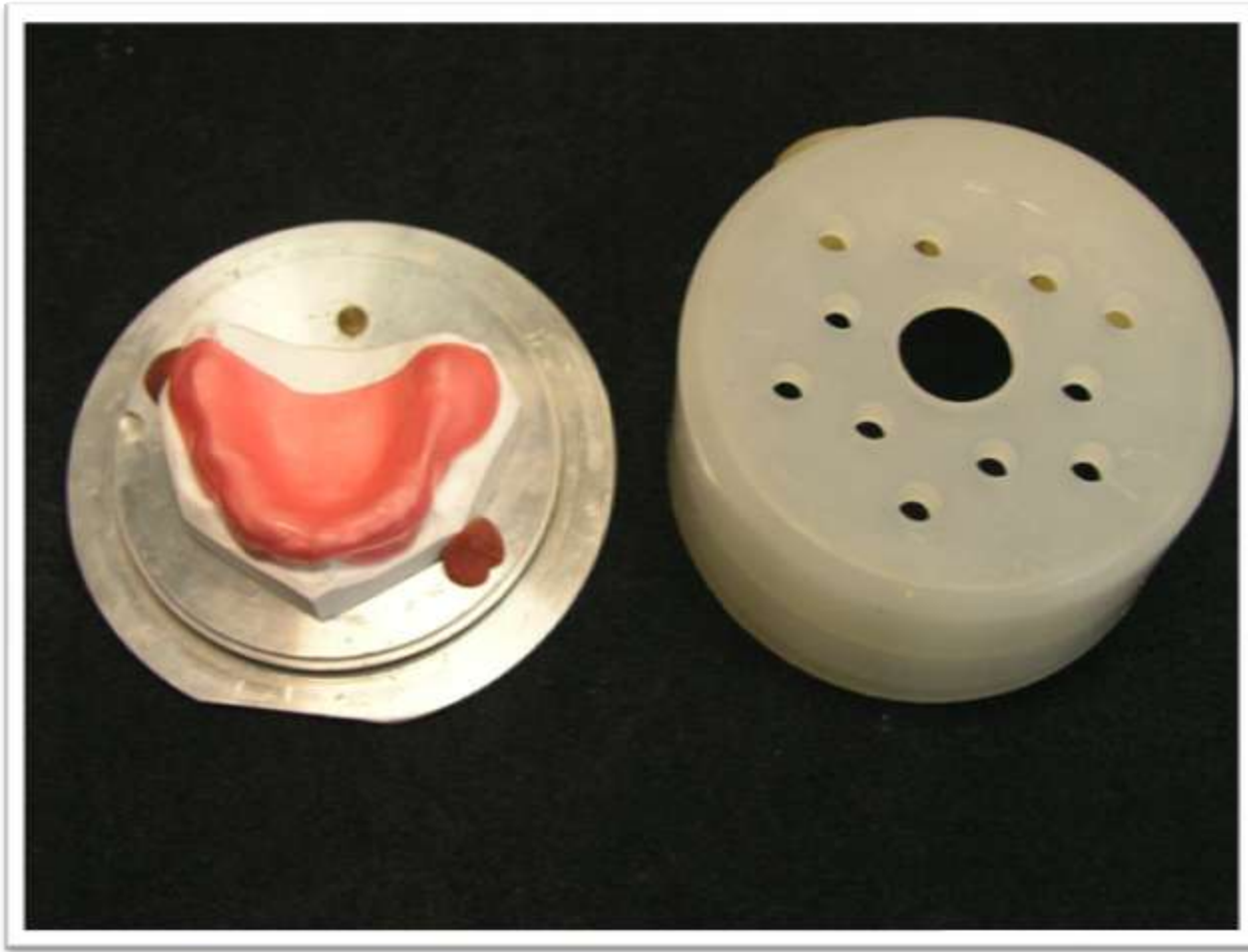

Type 2




Type 3

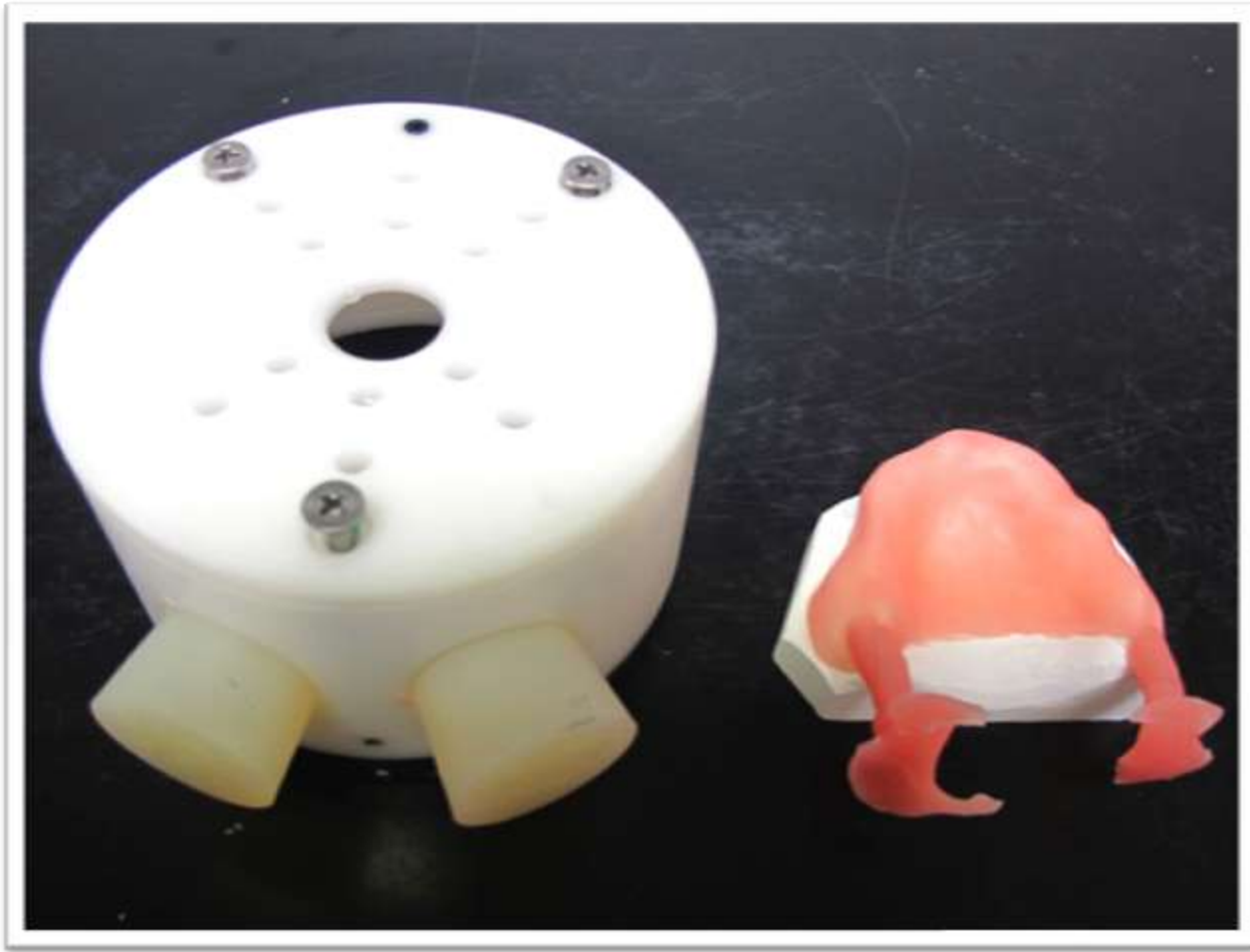


Type 3




Type 3

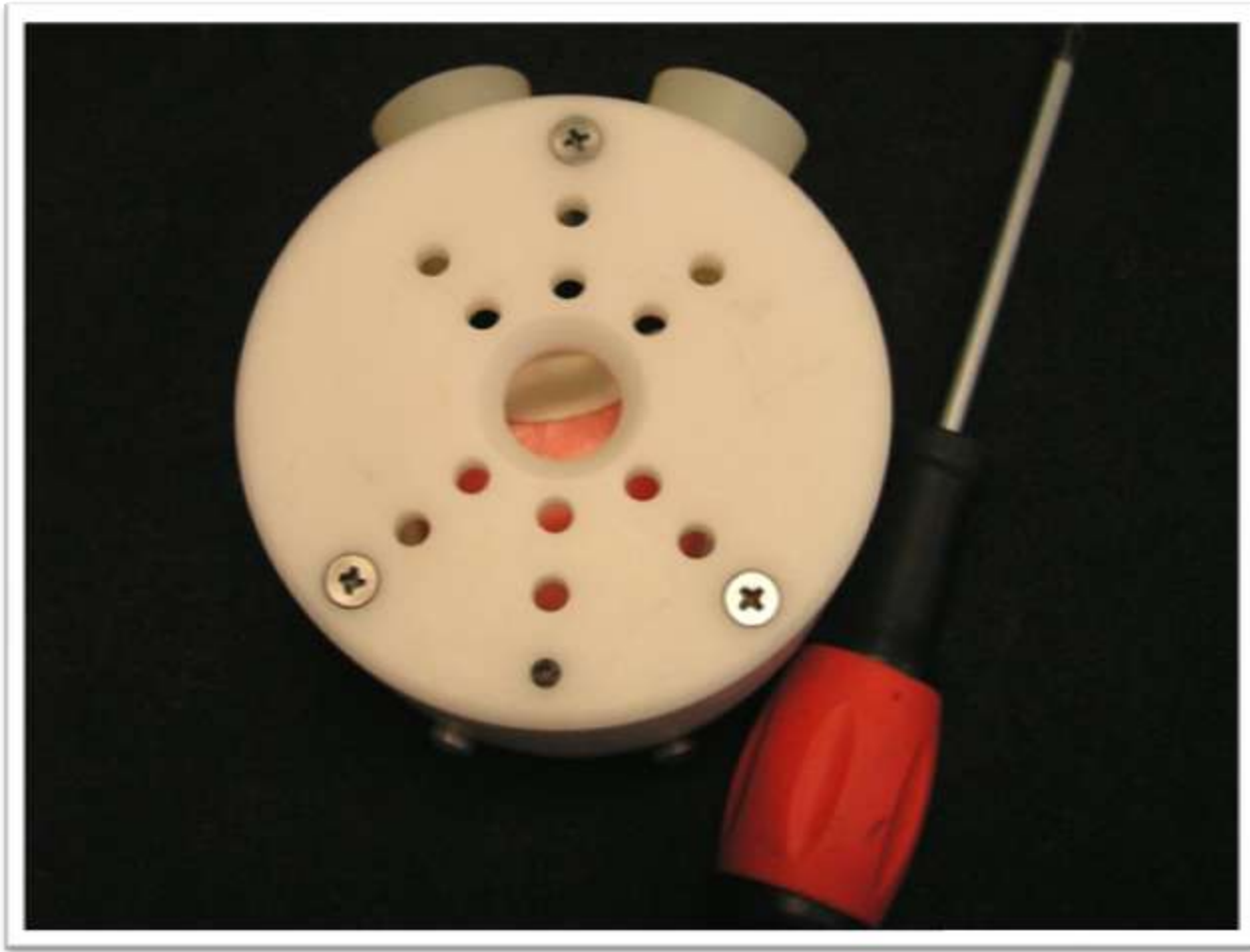

Type 4

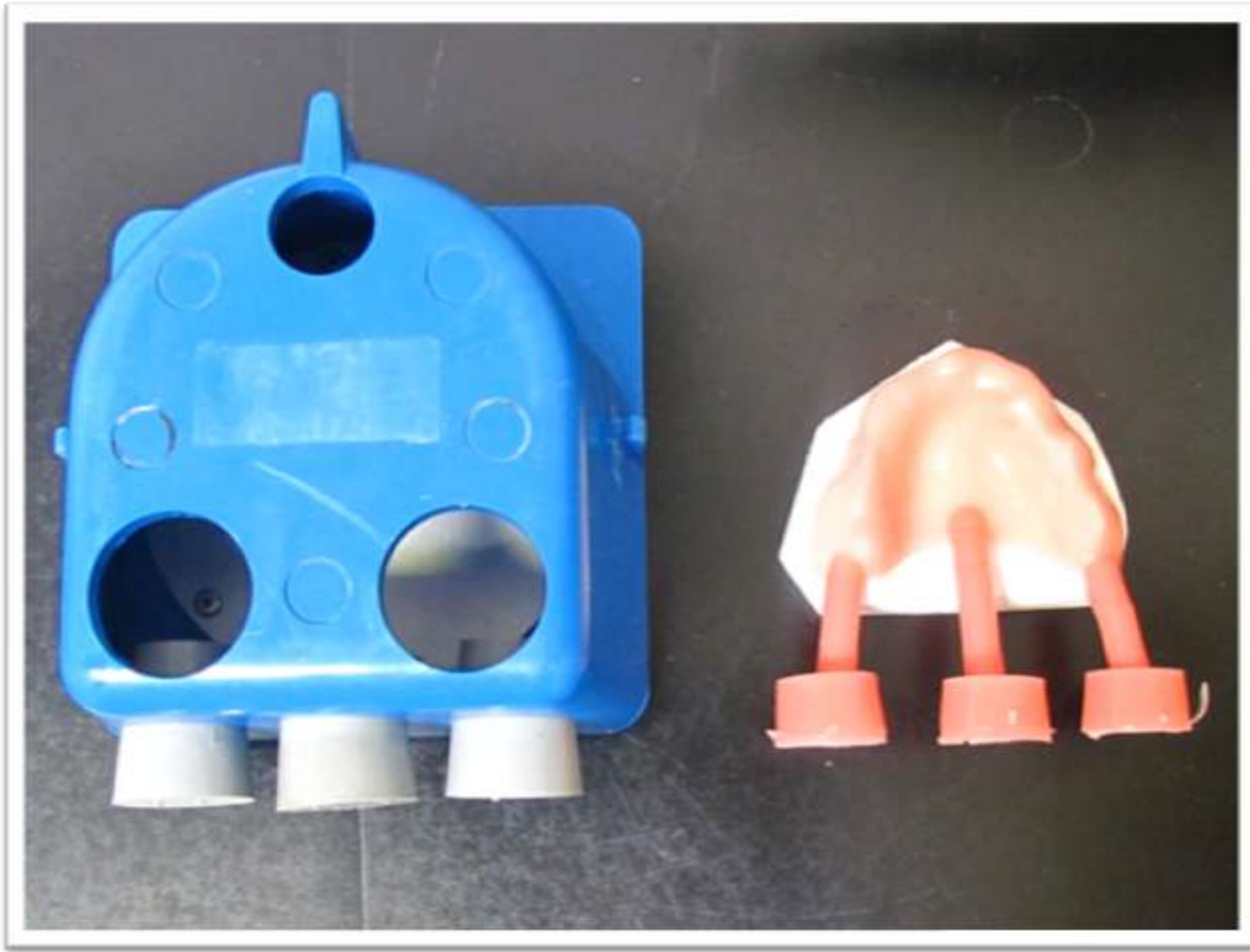

Type 4

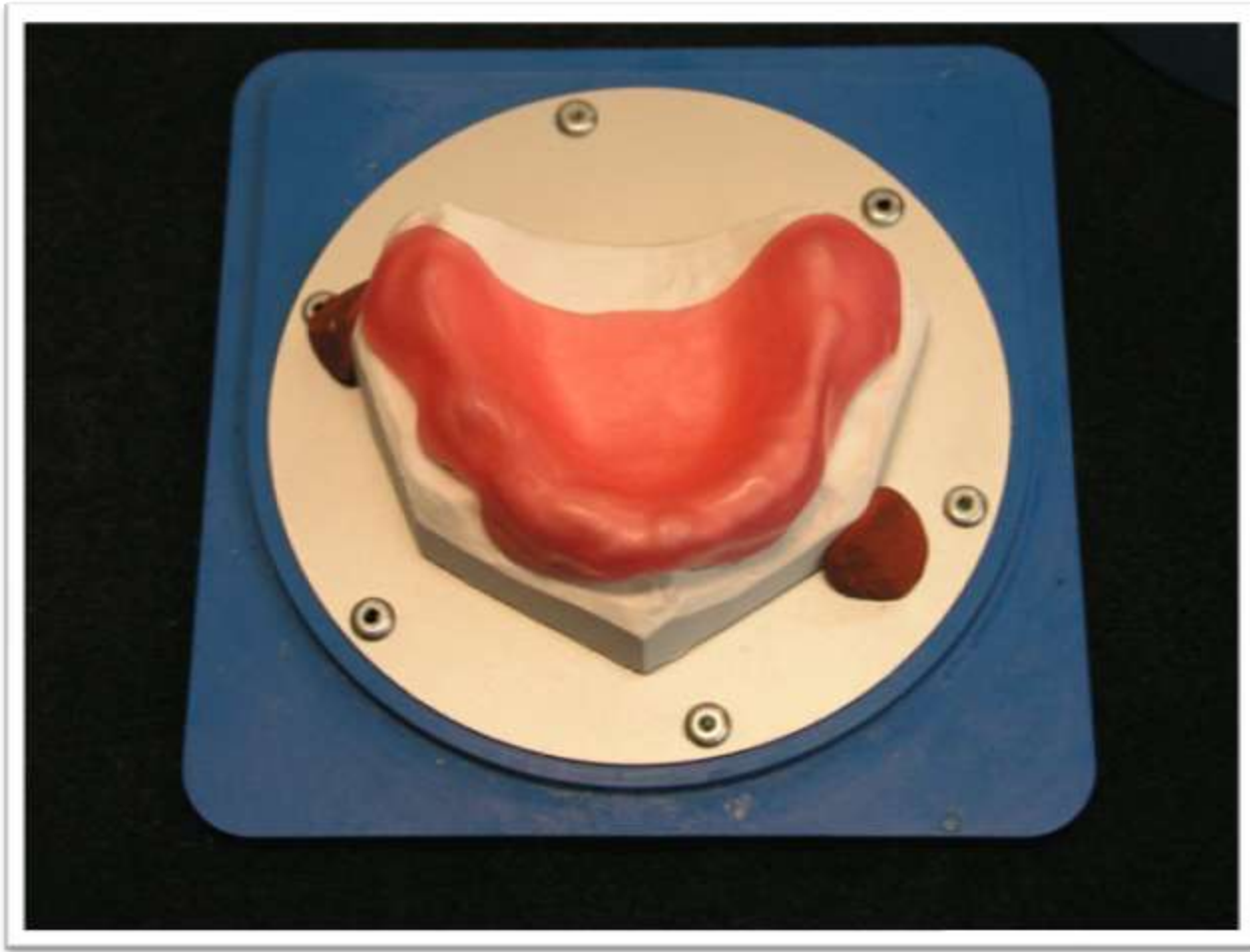

Type 4



Used quantity of duplicating materials

Flask	Duplicating material
MKM Transparent	250g
Type 1	290g
Type 2	300g
Type 3 (2 canals)	330g
Type 4	410g

One model was used for all measurements .

Number of castings per 6kg of used duplicating material

Flask	Castings
MKM Transparent	24x
Type 1	20x
Type 2	20x
Type 3 (2 canals)	18x
Type 4	14x

One model was used for all measurements.

Comparasion of used wax and resin

Flask	Resin	Wax
MKM Transparent	8,5g	4,3g
Type 1	8,6g	4,3g
Type 2	9,6g	4,2g
Type 3 (2 canals)	8,3g	4,1g
Type 4	11,0g	4,1g

One model was used for all measurements.

**High quality products for
your dental laboratory**



Contact

Haanova 31, 851 04 Bratislava

tel: +421 2 622 48 942

fax: +421 2 622 48 931

e-mail: info@mkmsystem.sk

www.mkmsystem.sk

